
# Real Time Measurement of Amplitude Change and Velocity Change

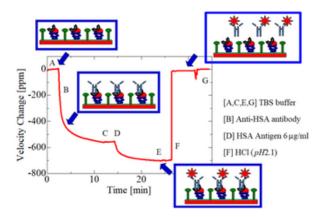


# **SAW Chip Specifications**

Features	Description
Size (LxWxH) / Weight	23 x 10 x 2 mm (Typ.) / 0.47 grm (Typ.)
Center Frequency	250 ± 2 MHz
Delay Time	1.2 μs ± 0.2 μs
Insertion Loss (bare chip)	24.2 dB - 32.5 dB
SAW Chip Substrate	Quartz
Sensor Surface	Gold
Pipetting Volume	5 μL (Typ.)
Sample Area Black Resin	Ероху
Sensor Channels	Two

#### **Notice for SAW chip**

- Use within 12 hours after opening a plastic package
- Do not touch the sensor surface
- Do not reuse SAW chips
- · Do not contaminate the electric terminal area
- Dispose of as industrial waste in accordance with the regulations



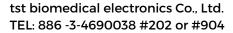







# tstbio SAW Biosensor

- Immunoassay analyzer
- Molecular internation for your research
- Development of your products



### **SAW Open Platform**



- Reader (with USB cable)
- Software on PC
- Protocol for immobilization
- SAW bare chip

#### **Features**

- 250 MHz SH-SAW measurement
- Small sample volume (3~5 μL)
- Manual pipetting
- Palm-sized Reader

#### **Notice for Platform Reader**

- Wait 10 minutes after reader turn on for warm up (check it at your test environment)
- · Place on a flat surface
- Do not contaminate the electric terminal area

## **Platform Reader Specifications**

Features	Description
Size (LxWxH)	107 x 73 x 30 mm (Typ.)
Weight	185 gram
Battery	Four AAA alkaline batteries
Signal Detection / Frequency	RF burst signal / 250MHz
Modes:	Modes:
1.Stand Alone Operation	1. Measure concentration
2. Real Time Measurement	2. Measure velocity (phase) and amplitude
QR Code Reader	Need to use for stand alone operation
Measurement Interval	0.125 seconds for 2CH
Number of Test Channels	Standard: two channels
Operating Conditions	Temperature: 25 ± 10°C,Humidity: 40 ± 30%
Software for PC	Platform Viewer Ver. 2.0.2.1
Software Version	Ver. 2.01.214C
SAW Chip	Bare chip (provided by tstbio) / Other (provided by tstbio)

#### **Gold Surface for Immobilization**

